

Optum

From the Office of Health Informatics **Oracle Health (Cerner) Millennium BCMA Self Administered Treatments**

September 19, 2025

Documenting the administration of patient self administered medication treatments within the Medication Administration Record (MAR) ensures a complete and accurate record. Timely and precise charting protects the patient from potentially dangerous errors, provides a clear communication tool for the entire care team, and establishes a legal record of care delivered.

Overview of Barcode Medication Administration

Barcode Medication Administration (BCMA) is a process designed to electronically confirm the 7 Patient Rights of Medication Administration: Right Patient, Right Medication, Right Dose, Right Time, Right Route, Right Documentation, and Right Reason. BCMA adds an additional layer of safety but DOES NOT replace the safe medication administration process. Hospital policies regarding Patient Identification (2-Identifier process) **MUST** be followed.


Self-Administered Treatment Administration Documentation Workflow

It may be appropriate for a patient to self-administer medication(s), such as an inhaler or nebulizer. When documenting this process, healthcare staff must follow established protocol to ensure accuracy and compliance. This includes observing the patient's self-administration and correctly entering the information into the Medication Administration Record (MAR).

- Prepare medication(s) to be administered. STEP 1:
 - If preparation is required, keep the barcode intact for scanning purposes.
- STEP 2: At the patient's bedside, verbally/visually verify the correct patient, following hospital policies.
- **STEP 3**: Confirm the patient information in the demographics banner bar.
- **STEP 4:** Open the **Medication Administration Wizard** from the barcode button on the toolbar.

IIIII Medication Administration

STEP 5: Scan the FIN barcode on the patient's wristband or face sheet per practice policy.

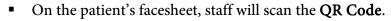



When scanning the QR Code, staff will need to scan in a downward fashion. This will avoid capturing the other barcodes on the sheet.

STEP 6: Click the **Refresh** button within the Wizard to ensure the most up-to-date medication order(s) are viewed.

Last Refresh at 13:52 EDT

- **STEP 7:** Scan the medication(s).
 - Complete any required or missing fields and manage alerts, as applicable.

• After scanning the patient wristband and medication(s), click in the **Result** column.



- The Charting For: window will open. In the
 *Performed by field, type Patient and
 PATIENT, SELF ADMIN will auto populate.
- Document any other required or pertinent information.
- Click **OK**.



Charting for: TESTING, NEILE

albuterol (Ventolin HFA)

2 Puffs, Aerosol, INHALATION, 05/08/20 15:00:00 EDT, 05/08/20 15:00:00 EDT, Treatment Location:
INLAND HOSPITAL
General Recommended Maximum Daily Dosing Guidelines: Nebulized Dose: Bronchospa

_ 🗆 X

STEP 8: Observe the patient self-administering the medication(s) to ensure appropriate administration.

STEP 9: Click the **Sign** button to document.