

The Clinical Leader Organizer (CLO) has icons that fire based on orders and nursing documentation within the EHR. These icons will assist leaders and nursing staff in identifying if documentation for certain measures has been completed as well as pinpointing certain risk factors.

**What do the Icons indicate and how do they populate on the Clinical Leader Organizer?**

- Click the cell to open the reading pane for more documentation or order details.
- Periodically refresh the CLO to get the most recent data. It does not automatically refresh.

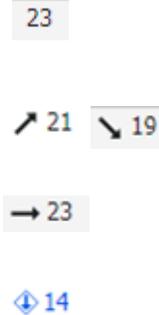
**NOTE:** These icons may not match the icons for the same thing in Capacity Management.

COLUMN	ICON OR INDICATOR	DESCRIPTION
Airway	 Tracheostomy Tube Activity:  Endotracheal Tube Activity:  Pharyngeal Airway Activity:	<p>This icon displays with the <b>type of airway</b> when the patient has an artificial airway documented.</p> <ul style="list-style-type: none"> <li>• Click the cell to view the activity, date, and time of last assessment, and who documented the assessment.</li> </ul>
Allergy	 Dust, No Known Medication Allergy  Darv... doxycy... am... clari... met... sim... p...	<p>This icon displays when the patient has documented <b>allergies</b>.</p> <ul style="list-style-type: none"> <li>• If the patient has multiple allergies, the allergies will be truncated.</li> <li>• Click the cell to view the complete list of allergens, reaction, and severity documentation.</li> </ul>
Allergy	 <b>No Allergies Recorded</b>	This icon displays when <b>No Allergy</b> documentation has been recorded.
Allergy	 No Known Medication Allergies	This icon displays when <b>No Known Medication Allergies</b> are documented.
Assessment	4-Eyed Skin Assessment: Yes Bedside Mobility Assessment: --  4-Eyed Skin Assessment: --  Bedside Mobility Assessment: Level	<p>Displays <b>4-Eyed Assessment</b> with a <b>Yes</b> if this assessment has been completed.</p> <p><b>Bedside Mobility Assessment</b> will display the <b>Level</b>.</p> <ul style="list-style-type: none"> <li>• -- will display if the assessment has not been completed.</li> </ul>

COLUMN	ICON OR INDICATOR	DESCRIPTION
Care Team	No Icon	<p>Click in the <b>Care Team</b> column to see who has established a relationship with the patient.</p> <ul style="list-style-type: none"> <li>• Click <b>Display Physician</b> to have the provider's name populate in the column for all patients on the list.</li> <li>• <b>Display Nursing</b> is available to populate the column once Display Physician is selected.</li> <li>• Nurse names will only appear under Assigned Nurses if the nurse is assigned using <b>Staff Assignment in Clairvia</b>.</li> <li>• Only one name will populate the column, a physician, or a nurse.</li> </ul> <p><b>NOTE:</b> The Nurse Tech name will display if the tech name comes in the alphabet before the nurse's name.</p>
Catheter	 Urethral Indwelling/Continuous 14 French	<p>This icon displays when the Urinary Catheter dynamic group has been created and Insert New Catheter or Daily Necessity Assessed has been documented in Activity.</p> <ul style="list-style-type: none"> <li>• Click the icon to see which Activity was documented.</li> </ul>
Catheter	 <b>Overdue</b> Urethral 3-Way 22 French	<p>This icon displays when no documentation in Activity has been documented in the last 72 hours.</p>
Central Line	 Antimicrobial coated catheter Triple Left, Upper Chest	<p>This icon displays when an active order for a Central Line is present.</p> <ul style="list-style-type: none"> <li>• It also displays when Activity and dressing change is documented.</li> </ul>
Central Line	 <b>Overdue</b> Peripherally inserted central catheter (PICC) Single Left, Bra	<p>This icon displays if Central Line dressing documentation has exceeded 7 days (168 hours).</p>
Discharge	 Discharge Ordered Order: -- Anticipated Discharge: JAN 06, 2025 00:00 Ordered On: --	<p>This icon displays when the patient has an Anticipated Day of Discharge documented but does not have a discharge order.</p>
Discharge	 <b>Discharge Overdue</b> Order: Discharge (DX) Anticipated Discharge: JAN 03, 2025 00:00 Ordered On: JAN 03, 2025 09:31	<p>This icon displays when the patient has a <b>Discharge order</b> placed.</p> <ul style="list-style-type: none"> <li>• Click the cell to view Order details.</li> </ul>

COLUMN	ICON OR INDICATOR	DESCRIPTION
Discharge	 Discharge Completed Order: -- Completed On: JAN 01, 2025 10:50	This icon displays once the <b>discharge</b> has been completed. <ul style="list-style-type: none"> <li>If using a location list, the patient will fall off the list once discharged and the CLO is refreshed, so this may not be seen unless a list for discharged patients is being used.</li> </ul>
Elopement	 Elopement Risk	This icon will only display for <b>Acadia</b> if the nurse has documented in <b>Environmental Safety Management &gt; BH Precautions Observed</b> both <b>Initiate</b> and <b>Elopement risk order</b> .
Falls Risk		This indicator displays the most recently documented <b>fall risk score</b> . <ul style="list-style-type: none"> <li>If only one assessment has been documented, just a number will display.</li> <li>A trending up or down arrow will indicate the score trend.</li> <li>If no change from the previous value, an arrow will display horizontally.</li> <li>Click the cell to see the last assessment date and who completed it.</li> </ul>
Generic High Risks Column	 2 High Risk Documented	This icon displays with a number indicating how many high-risk orders have been placed for the patient. <ul style="list-style-type: none"> <li><b>Precautions (Nsg)</b>- any of the options in the order, <b>Continuous Observation</b>, or <b>Room Search (Nsg)</b> will cause the icon to populate.</li> <li><b>Suicide Precautions</b></li> <li><b>Elopement Precautions</b> (Acadia only)</li> <li>Click the cell to see the type of precaution &amp; the order date.</li> </ul>
Location	No Icon	The location column displays the <b>nursing unit name</b> , <b>room number</b> , and <b>bed placement</b> .
Medication History		This icon displays when the <b>Medication History</b> has not been started.

COLUMN	ICON OR INDICATOR	DESCRIPTION
Medication History		This icon displays when the <b>Medication History</b> has been completed. <ul style="list-style-type: none"> <li>Click the icon to see when it was completed and by whom.</li> </ul>
Observation	OBS	This indicator displays when the patients <b>encounter status</b> is <b>Observation</b> . No data displays if icon is clicked.
Oxygen Therapy	 	This indicator will display if <b>Initiate</b> is documented in <b>Oxygen Activity</b> . <ul style="list-style-type: none"> <li>If oxygen is documented and no flow rate is documented, - - will display for Rate.</li> <li>Click the cell to see the FIO2, when and by whom this information was documented.</li> </ul>
Pain	   	This indicator displays the most recently documented <b>pain score</b> . <ul style="list-style-type: none"> <li>If only one assessment has been documented, just a number will display.</li> <li>A trending up or down arrow will indicate the score trend.</li> <li>If no change from the previous value, an arrow will display horizontally.</li> <li>Click the cell to see the last three pain scores documented.</li> </ul>
Patient	No Icon	The Patient column displays the <b>patient's name, age, and gender</b> . <ul style="list-style-type: none"> <li>Click the patient's name to <b>open</b> the patient <b>chart</b>.</li> <li>Click the cell to view additional data such as <b>DOB, MRN, and FIN</b>.</li> </ul>
Restraints	 Behavioral Restraint  Non Behavioral Restraint	This icon displays when an order for Behavioral or Non-Behavioral Restraints has been placed. <ul style="list-style-type: none"> <li>Click the cell to see the order details.</li> </ul>
Resuscitation	No Icon	No icon displays when the patient has a code status of <b>Full Resuscitation</b> . Click the cell for order details.

COLUMN	ICON OR INDICATOR	DESCRIPTION
Resuscitation		This icon displays when the patient has a Code Status and Treatment Restriction order for <b>DNAR/DNI</b> and <b>DNAR/Comfort Care</b> . Click the cell for order details.
Resuscitation		This icon displays when the patient has a Code Status and Treatment Restriction order for <b>DNAR/May Intubate</b> .
Skin Integrity		<p>This indicator displays the most recently documented <b>Braden Score</b>.</p> <ul style="list-style-type: none"> <li>• If only one assessment has been documented, just a number will display.</li> <li>• A trending up or down arrow will indicate the score trend.</li> <li>• If no change from the previous value, an arrow will display horizontally.</li> <li>• If the score is <b>below 18</b>, an <b>Out of Range Low</b> diamond will display in blue to the left of the number.</li> <li>• Click the cell to view the last 3 documented scores.           <ul style="list-style-type: none"> <li>○ Click <b>View Details</b> to view documentation results and by whom.</li> </ul> </li> </ul>
Suicide		This icon displays <b>if Suicidal Ideations</b> are documented, and the patient does not have a Suicide Precautions order.
Telemetry		<p>This icon indicates the patient has a <b>Telemetry Monitoring order</b>.</p> <ul style="list-style-type: none"> <li>• Click the cell to view the order details.</li> </ul>
Ventilator		<p>This icon displays when there is a <b>Vent Settings order</b>.</p> <ul style="list-style-type: none"> <li>• Click the icon to view the order details.</li> </ul>