


Accurate documentation of vital signs from FetaLink into PowerChart ensures a complete record for the patient. Frequent documentation during procedures will be addressed in this flyer.

## Frequent Vitals Signs

Through BMDI, the frequency of vital signs can be set to as frequently as needed. The frequency for vital signs in iView needs to be adjusted to match those.

### ➤ **Setting the timeframe in iView**

- Within iView, near the date and time, right-click to select Q1 min.
- Select **Clinical Range**. The time range will default to 8 hours and select **OK**.

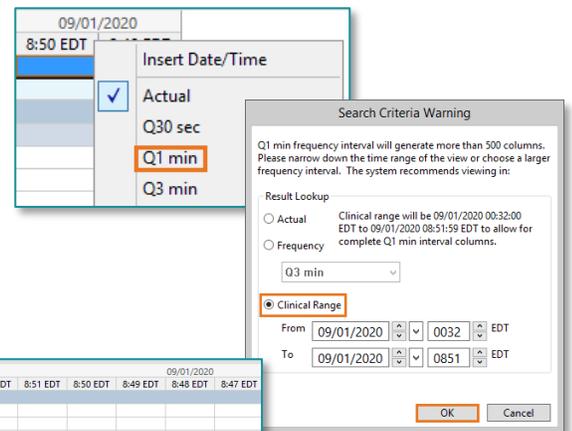

### ➤ **Verification of Vital Signs**

- Vital Signs will pull in automatically in an unvalidated state (pink text) from FetaLink to iView.
- Double click on the time column of the vital signs after verifying them to change them to a validated state (the text will change purple.)
- Multiple columns may be validated before signing them (the text will change to black).
- Click the green checkmark to sign the results.

Vital Signs	8:54 EDT	8:53 EDT	8:52 EDT	8:51 EDT	8:50 EDT	8:49 EDT	8:48 EDT
Temperature							
Temperature Fahrenheit							
Temperature Method							
Cuff SBP/DBP					116/79		118/76
Cuff MAP Measured					110		112
Cuff MAP Estimated					110		
Cuff Location							
Cuff Method							
Apical Heart Rate					82		79
Pulse Rate							
Pulse Regular							
Heart Rate-Monitored							
Cardiac Rhythm							
Respiratory Rate							
Oxygen Activity							
Oxygen Flow Rate							
Oxygen Delivery							
Oxygen Saturation					100	99	100
Oxygen Sat Activity							
Supine SBP/DBP							
Supine Pulse							
Sitting SBP/DBP							
Sitting Pulse							
Standing SBP/DBP							
Standing Pulse							




Vital Signs	8:54 EDT	8:53 EDT	8:52 EDT	8:51 EDT	8:50 EDT	8:49 EDT	8:48 EDT
Temperature							
Temperature Fahrenheit							
Temperature Method							
Cuff SBP/DBP					116/79		118/76
Cuff MAP Measured					110		112
Cuff MAP Estimated					91		90
Cuff Location							
Cuff Method							
Apical Heart Rate					82		79
Pulse Rate							
Pulse Regular							
Heart Rate-Monitored							
Cardiac Rhythm							
Respiratory Rate							
Oxygen Activity							
Oxygen Flow Rate							
Oxygen Delivery							
Oxygen Saturation					100	99	99
Oxygen Sat Activity							
Supine SBP/DBP							
Supine Pulse							
Sitting SBP/DBP							
Sitting Pulse							
Standing SBP/DBP							
Standing Pulse							



**NOTE:** Only purple verified vitals will be signed. Pink vital signs will not be recorded in the chart.

➤ **Changing the time frame when frequent vital signs are not needed.**

- From iView, right-click in the time column and select **Actual time**.
- Review the time of recorded vitals from the **FetaLink annotations summary**.
  - You can Select All or only select Vital Signs for a filtered view.
- Right click to set the appropriate time in iView closest to the nearest minute of the time shown in FetaLink.
- Double click on the time to pull in the vitals from FetaLink.
- Verify the vitals and sign.

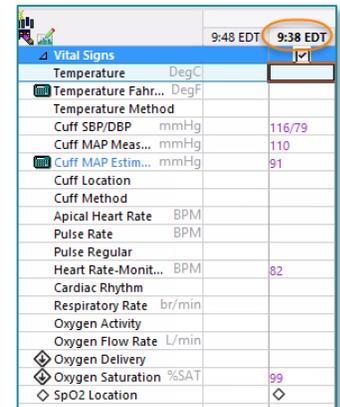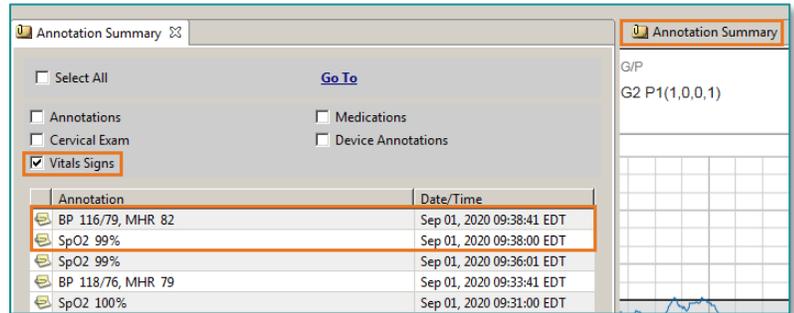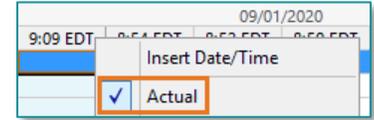

➤ **Changing the clinical timeframe**

If the timeframe needed to pull vitals into is not in view, the clinical range can be changed.

- From the timeframe toolbar, right click and select a timeframe option or click Other.
- When selecting Other, a specific timeframe may be selected.

