




Cerner Millennium provides functionality that allows end-users to view radiographic images within the electronic medical record (EMR). See below for details.

Viewing Radiographic Images in Millennium



- Open the desired imaging report in one of the following areas:
 - Summary/Workflow MPages
 - Notes
 - Results Review
 - Documentation
- Click the View Image  icon in the toolbar at the top of the report.

NOTE: The View Image icon is dithered out initially () and is only available when you have opened a report. ()



Tools Available in the Millennium Radiology Viewer

- To access available tools, you will need to hover over the left hand side of the image workspace (indicated below by the 1:).
- Adjust Brightness/Contrast (Window Level)
 - This is the default tool.
 - Invert
 - Original Window Level

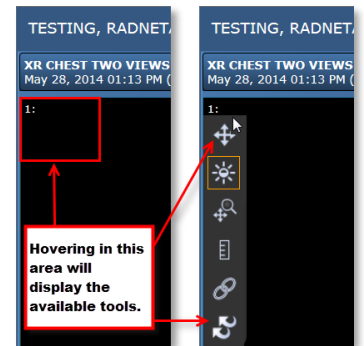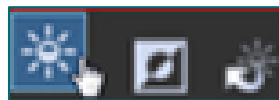


NOTE: To adjust the brightness of the image, click and drag (left to right; top to bottom) on the image itself. This function is different than IMPAX!


- Image Navigate /Flip/Rotate
 - Rotate right
 - Flip Horizontal



- Zoom/Pan - Use the scroll wheel to zoom in or out from the current position of the mouse pointer.
 - Zoom
 - Magnify
 - Zoom In
 - Zoom Out
 - Magnify Pan Reset

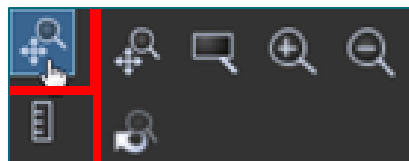

Toggle Manual Link

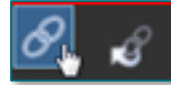

NOTE: You will need to link both images.


➤ Reset Manual Series Linking

- Reset (this icon will reset the entire image back to its original view).



➤ Options Menu

- Click the Options  icon to access Study Comments and IMPAX Reports:

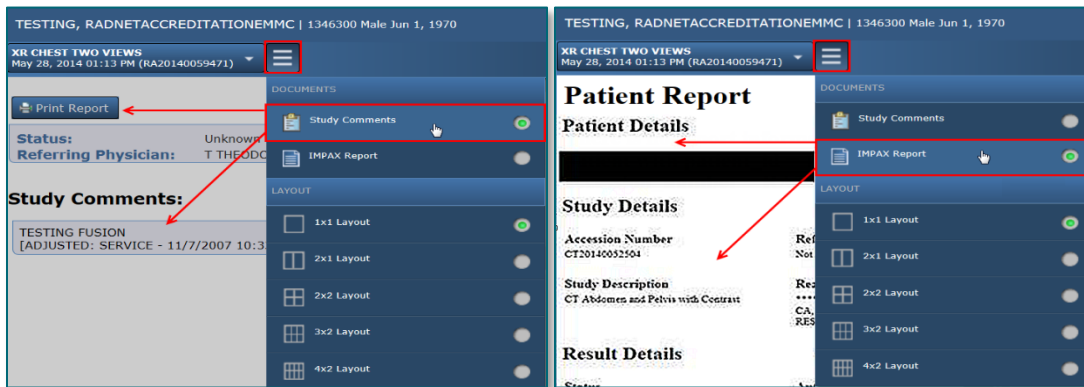

- The Study Comments option is only viewable if there are actual comments associated with the image. If there are no comments, Study Comments will not display.

➤ Located to the right are additional **quick link icons** (e.g. Reset, Navigate Previous, Navigate Next, Display Series, Compare Studies).

